

4 Systems to Choose From

- A complete selection of large-format submarine gel systems
- Each system is 200 mm wide and range in length from 200 mm to 400 mm
- Designed for maximum flexibility, combs interchange between models and shorter gel trays may be used in longer format chambers


Simple Gel Casting

- Gasketed gel tray end gates ensure ease of handling (Model CSSU2027 uses tape)

- Press the end gates on each gel tray cast the gel and pull them off
- Each tray is UV transparent

Easy Sample Loading

- Each tray accommodates multiple comb positions, allowing you the option of running two or more rows of samples on every gel
- Each tray is equipped with red visualisation stripes to aid sample loading

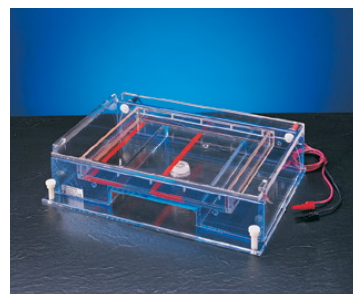

Thermo Electron's Classic Series Maxigel Gel Systems come in four convenient sizes 200 x 200 mm, 200 x 250 mm, 200 x 270 mm and 200 x 400 mm.


Ordering Information

Cat No	Description
CSSU2020-LVD	Includes: main chamber, 200 x 200 mm gel tray, 2 x gel tray end gates, 1 x 24/36 well 1.5 mm double-sided comb, safety cover with leads and instructions
CSSU2020-115	Includes: main chamber, 200 x 200 mm gel tray, 2 x gel tray end gates, 1 x 24/36 well 1.5 mm double-sided comb, safety cover with leads and instructions
CSSU2025-LVD	Includes: main chamber, 200 x 250 mm gel tray, 2 x gel tray end gates, 1 x 24/36 well 1.5 mm double-sided comb, safety cover with leads and instructions
CSSU2025-115	Includes: main chamber, 200 x 250 mm gel tray, 2 x gel tray end gates, 1 x 24/36 well 1.5 mm double-sided comb, safety cover with leads and instructions
CSSU2040-LVD	Includes: main chamber, 200 x 400 mm gel tray, 2 x gel tray end gates, 1 x 24/36 well 1.5 mm double-sided comb, safety cover with leads and instructions
CSSU2040-115	Includes: main chamber, 200 x 400 mm gel tray, 2 x gel tray end gates, 1 x 24/36 well 1.5 mm double-sided comb, safety cover with leads and instructions
CSSU2027-LVD	Includes: main chamber, 200 x 270 mm gel tray, 1 x 24/36 well 1.5 mm double-sided comb, safety cover with leads and instructions
CSSU2027-115	Includes: main chamber, 200 x 270 mm gel tray, 1 x 24/36 well 1.5 mm double-sided comb, safety cover with leads and instructions

Classic Series 96 Well Compatible Gel System

Product Features

- Developed specifically for multiple sample/high throughput analysis accommodating multichannel pipettes and 96 well plates
- Two gel tray lengths 100 mm and 150 mm
- Both are provided in this 250 mm wide chamber that utilises 25 well and 50 well combs
- Each comb is designed to accept samples directly from a multichannel pipette.
- The 25 well comb allows for 24 sample wells with one reference well
- The 50 well comb allows for 48 sample wells and 2 reference wells
- Each gel tray is designed to accept two combs for maximum sample throughput
- Gel casting is facilitated by the use of gasketed gel trays that press into a special casting stand
- Each gel tray is UV transparent and has two red stripes to aid sample loading

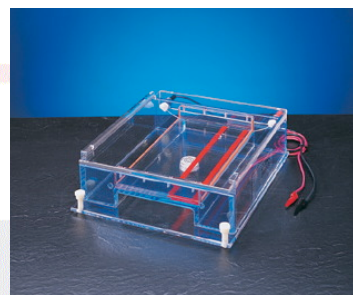

Ordering Information

Cat No	Description
CSSU2515-LVD	Includes: main chamber, 100 mm gel tray, 150 mm gel tray, 1 casting stand 1 x 25/50 well dual sided 1.5 mm comb, safety cover with leads and instructions
CSSU2515-115	Includes: main chamber, 100 mm gel tray, 150 mm gel tray, 1 casting stand 1 x 25/50 well dual sided 1.5 mm comb, safety cover with leads and instructions
CSSU2515-2515M	25 & 50 well double sided Microtiter format comb, 1.5 mm
CSSU2515-25MA	25 well Microtitre comb, 1.5 mm
CSSU2515-25MB	25 well Microtitre comb, 1.0 mm
CSSU2515-50MA	50 well Microtitre comb, 1.5 mm
CSSU2515-50MB	50 well Microtitre comb, 1.0 mm
CSST2510	25 x 100 mm gasketed gel tray with 2 red stripes
CSST2515	25 x 150 mm gasketed gel tray with 2 red stripes
CSST2515GK	250 mm gasket kit (pair)
CSSC2515	Casting stand (trays not included)