

Applications

For surgical sterilisation of diverse instruments
surgical odontological, etc

Adjustable temperatures from 60°C up to 250°C. stability +/-6 °C

Features

Heating by shielded elements in the base which provide a rapid temperature rise
Flap door.
Inner chamber in AISI 304 stainless steel with a heater cover, three shelf runners and two perforated shelves with 10 mm high.
Epoxy-coated outer casing

Safety

Safety thermostat incorporated EN.61010 Standard.

Control Panel

Hydraulic thermostat for temperature control.
Locking device for thermostat knob.
Timer 0 to 120 min. with automatic off.
Analogue temperature reading thermometer.



Part No.	Capacity litres	H x W x D (int. dims) cm	H x W x D (ext. dims) cm	Power W	Weight Kgs
C/O/010002	19	25 x 32 x 23	37 x 54 x 34	770	20

Glass Drying Oven Dryglass

Features

Hydraulic thermostat for temperature control.
Air circulation by turbo fan.
Inner chamber made of AISI 304 stainless steel with shelf runner devices.
Removable tempered glass sliding doors.
Ventilation holes for steam. Epoxy coated external case.

Safety

Over temperature safety cut out incorporated EN.61012 Standard.
Adjustable over temperature safety thermostat can be added as an accessory to comply with the DIN 12.880.2 Standard.


Standard equipment

2 shelves and 4 guides.

Control Panel

Dual heating power selector switch.
Mains indicator lamp.
Hydraulic thermostat for temperature control.
Locking device for thermostat knob.
Heater operation indicator lamp.
Analogue temperature reading thermometer.

Adjustable temperatures from 40°C up to 70°C. stability +/-6 °C

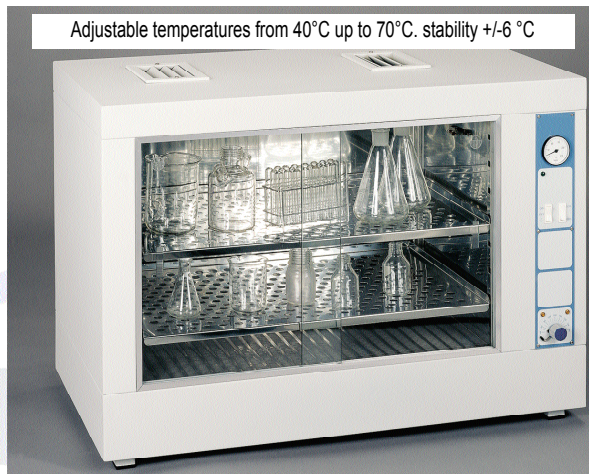

Part No.	Capacity litres	H x W x D (int. dims) cm	H x W x D (ext. dims) cm	Shelf Positions	Power W	Weight Kgs
C/O/010003	126	45 x 70 x 40	66 x 94 x 54	8	3000	58

Accessories (Accessories must be factory installed)

C/O/010100 - Safety thermostat which disconnects the heating in case of oven's own thermostat failing, comes with manual reset (conforms to DIN 12.880 class 2)
C/O/010101 - Timer switch 0-120 minutes.
C/O/010102 - Timer switch 0-12 hours.
C/O/010103 - 24 hour programmer with continuous on/off cycling up to every 15 minutes.

Spares

C/O/010200 - Shelf guides x 2.
C/O/010201 - Shelf.
Each shelf requires 2 guides.